

# Human MICA ELISA Instructions

**Cat: EH0044**

## Content

	CAT	Volume
① CP (Coated Plate)	EH0044CP	96 well
② S (Standard)	EH0044S,S1-S7,S0	9 vial
③ AB (Assay Buffer 1×)	EAB01	12 ml/bottle
④ SD (Sample Diluent)	ESD03	12 ml/bottle
⑤ DA-H (Detect Antibody-HRP)	EH0044DA-H	6 ml/bottle
⑥ TS (TMB Substrate)	ETS01	12 ml/bottle
⑦ SS (Stop Solution)	ESS01	12 ml/bottle
⑧ WB (Wash Buffer 10×)	EWB01	50 ml/bottle
⑨ SF (Sealer Film)	ESF01	6 pieces

**NOTE:** After the kit is opened, the stabilization period of each content is 30 days, so please use the kit within 30 days after opening.

## REAGENT PREPARATION

### Washing Buffer (1×) Preparation


Pour entire contents (50 ml) of the **Washing Buffer Concentrate** (10×) into a clean 500 ml graduated cylinder. Bring to final volume of 500 ml with glass-distilled or deionized water. Transfer to a clean wash bottle and store at 2 to 25°C.

### Standard Curve Preparation:

S1 to S7 and S0 is ready to use for serum and plasma.

Other sample type, prepare the standard curve with whatever buffer (SPB, Sample Prepared Buffer) is used to prepare the sample, such as cell culture supernatant, tissue grinding liquid, cell lysate, etc. Urine sample use AB (Assay Buffer) prepare standard curve.

The human MICA Standard EH0044S 100000 pg/ml 40  $\mu$ l + 360  $\mu$ l SPB serves as the high standard (10000 pg/ml). Pipette 200  $\mu$ l of SPB into each tube. Use the high standard to produce a 1:1 dilution series. Mix each tube thoroughly before the next transfer. SPB serves as the zero standard (0 pg/ml).

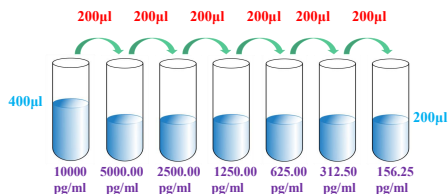

## ASSAY PROCEDURE

**Bring all reagents and samples to room temperature before use.**

- ① Prepare all reagents and working standards as directed in the previous sections.
- ② Remove excess **CP** (Coated Plate) strips from the plate frame, return them to the foil pouch and reseal.

- ③ Add 50  $\mu$ l of **AB** (Assay Buffer) to each well.
- ④ Add 50  $\mu$ l of **Standard or sample** per well. Ensure reagent addition is uninterrupted and completed within 15 minutes.

- ⑤ Add 50  $\mu$ l of **DA-H** (Detect Antibody-HRP) to each well.


- ⑥ Cover with an **SF** (Sealer Film). Incubate at room temperature (18 to 25°C) for 1 hour on a microplate **shaker** set at 500 rpm.

- ⑦ Aspirate each well and **wash**, repeating the process four times. Wash by filling each well with **WB** (Washing Buffer 300  $\mu$ l). Complete removal of liquid at each step is essential to good performance. After the last wash, remove any remaining **WB** (Washing Buffer) by aspirating or decanting. Invert the plate and **blot** it against clean paper towels.

- ⑧ Add 100  $\mu$ l of **TS** (TMB Substrate) to each well. Incubate for 5-30 minutes at room temperature.

- ⑨ Add 100  $\mu$ l of **SS** (Stop Solution) to each well.
- ⑩ Determine the optical density within 30 minutes, using microplate reader set to 450 nm corrected with 570 nm or 630 nm.

## TYPICAL DATA

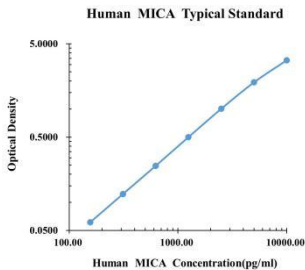

pg/ml	O.D.		Average	Corrected
0.00	0.0202	0.0214	0.0208	
156.25	0.0808	0.0828	0.0818	0.0610
312.50	0.1430	0.1440	0.1435	0.1227
625.00	0.2682	0.2661	0.2672	0.2464
1250.00	0.5125	0.5294	0.5210	0.5002
2500.00	1.0160	1.0370	1.0265	1.0057
5000.00	1.9670	1.9410	1.9540	1.9332
10000.00	3.2950	3.4180	3.3565	3.3357

## SENSITIVITY

The minimum detectable dose (MDD) of human MICA is typically less than 8.30 pg/ml.

The MDD was determined by adding two standard deviations to the mean optical density value of ten zero standard replicates and calculating the corresponding concentration.

## PRECISION

**Intra-assay Precision (Precision within an assay)** Three samples of known concentration were tested twenty times on one plate to assess intra-assay precision.

**Inter-assay Precision (Precision between assays)**

	Intra-assay Precision			Inter-assay Precision		
	S1	S2	S3	S1	S2	S3
Sample Number	22	22	22	6	6	6
Average (pg/ml)	389.9	1231.6	3586.9	373.1	85.4	3682.9
Standard Deviation	21.2	45.6	111.9	4.6	41.5	85.4
Coefficient Of Variation (%)	5.4	3.7	3.1	2.5	3.1	2.3

## RECOVERY

The spike recovery was evaluated by spiking 3 levels of human MICA into health human serum sample. The un-spiked serum was used as blank in this experiment.

The recovery ranged from 94% to 112% with an overall mean recovery of 106%.

## LINEARITY

To assess the linearity of the assay, five samples were spiked with high concentration of MICA in human serum and diluted with Sample Diluent to produce samples with values within the dynamic range of the assay.

The linearity ranged from 97% to 109% with an overall mean recovery of 103%.

## SAMPLE VALUES

Serum/Plasma – Thirty samples from apparently healthy volunteers were evaluated for the presence of MICA in this assay. No medical histories were available for the donors.

Sample Matrix	Sample Evaluated	Range (pg/ml)	Detectable %	Mean of Detectable (pg/ml)
Serum	30	469.25-3717.25	100	1351.05

n.d. = non-detectable. Samples measured below the sensitivity are considered to be non-detectable.