

# Mammal INHBE ELISA Instructions

## Cat:EHY0245

### Content

	CAT	Volume
① CP (Coated Plate)	EHY0245CP	96 well
② S (Standard)	EHY0245S1-S7,S0	8 vial
③ DA-H (Detect Antibody-HRP)	EHY0245DA-H	6 ml/bottle
④ AB (Assay Buffer 1×)	EAB01	12 ml/bottle
⑤ TS (TMB Substrate)	ETS01	12 ml/bottle
⑥ SS (Stop Solution)	ESS01	12 ml/bottle
⑦ WB (Wash Buffer 10×)	EWB01	50 ml/bottle
⑧ SF (Sealer Film)	ESF01	6 pieces

**NOTE:** After the kit is opened, the stabilization period of each content is 30 days, so please use the kit within 30 days after opening.


### REAGENT PREPARATION

#### ■ Washing Buffer (1×) Preparation

Pour entire contents (50 ml) of the **Washing Buffer Concentrate** (10×) into a clean 500 ml graduated cylinder. Bring to final volume of 500 ml with glass-distilled or deionized water. Transfer to a clean wash bottle and store at 2 to 25°C.

#### ■ Standard Curve Preparation:

S1 to S7 and S0 is ready to use for serum and plasma.

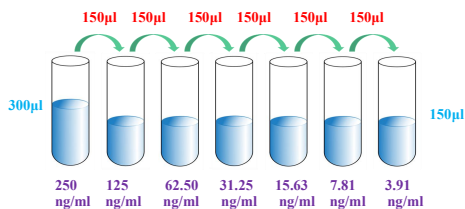


### ASSAY PROCEDURE

Bring all reagents and samples to room temperature before use.

- ① Prepare all reagents and working standards as directed in the previous sections.
- ② Remove excess **CP** (Coated Plate) strips from the plate frame, return them to the foil pouch and reseal.
- ③ Add 50 µl of **AB** (Assay Buffer) to each well.
- ④ Add 10 µl of **Standard or sample** per well. Ensure reagent addition is uninterrupted and completed within 15 minutes.
- ⑤ Add 50 µl of **DA-H** (Detect Antibody-HRP) to each well.
- ⑥ Cover with an **SF** (Sealer Film). Incubate at room temperature (18 to 25°C) for 60 min on a microplate shaker set at 500 rpm.
- ⑦ Aspirate each well and **wash**, repeating the process four times. Wash by filling each well with WB (Washing Buffer 300 µl). Complete removal of liquid at each step is essential to good performance. After the last wash, remove any remaining WB (Washing Buffer) by aspirating or decanting. Invert the plate and **blot** it against clean paper towels.
- ⑧ Add 100 µl of **TS** (TMB Substrate) to each well. Incubate for 5-30 minutes at room temperature.
- ⑨ Add 100 µl of **SS** (Stop Solution) to each well.
- ⑩ Determine the optical density within 30 minutes, using microplate **reader** set to 450 nm corrected with 570 nm or 630 nm.

## TYPICAL DATA

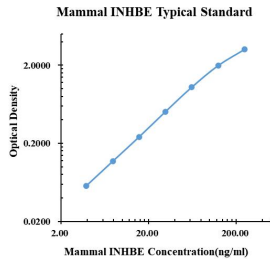

ng/ml	O.D.	Average	Corrected
0.00	0.0073	0.0065	0.0069
3.91	0.0656	0.0629	0.0643
7.81	0.1209	0.1289	0.1249
15.63	0.2501	0.2472	0.2487
31.25	0.5334	0.4981	0.5158
62.50	1.0980	1.0110	1.0545
125.00	2.0750	1.8970	1.9860
250.00	3.2990	3.0870	3.1930

## SENSITIVITY

The minimum detectable dose (MDD) of mammal INHBE is typically less than 0.06 ng/ml (10 µl of sample volume).

The MDD was determined by adding two standard deviations to the mean optical density value of ten zero standard replicates and calculating the corresponding concentration.

## PRECISION

**Intra-assay Precision (Precision within an assay)** Three samples of known concentration were tested twenty times on one plate to assess intra-assay precision.

**Inter-assay Precision (Precision between assays)**

	Intra-assay Precision			Inter-assay Precision		
	S1	S2	S3	S1	S2	S3
Sample Number	22	22	22	6	6	6
Average (ng/ml)	3.4	15.1	45.8	4.2	20.0	53.8
Standard Deviation	0.1	0.6	2.3	0.2	0.6	1.3
Coefficient of Variation (%)	2.6	3.7	4.9	4.7	3.2	2.4

## RECOVERY

The spike recovery was evaluated by spiking 3 levels of mammal INHBE into health mammal serum sample. The un-spiked serum was used as blank in this experiment.

The recovery ranged from 80% to 119% with an overall mean recovery of 96%.

## LINEARITY

To assess the linearity of the assay, five samples were spiked with high concentration of INHBE in mammal serum and diluted with Sample Diluent to produce samples with values within the dynamic range of the assay.

The linearity ranged from 81% to 117% with an overall mean recovery of 95%.

## SAMPLE VALUES

Serum/Plasma – Thirty samples from apparently healthy volunteers, thirty samples from apparently healthy mice, thirty samples from apparently healthy rats and two samples from healthy monkeys were evaluated for the presence of mammal INHBE in this assay. No medical histories were available for the donors.

Sample Matrix	Sample Evaluated	Range (ng/ml)	Detectable %	Mean of Detectable (ng/ml)
Human Serum	30	3.22-282.5 4	100	30.90
Mouse Serum	30	6.54-28.08	100	13.48
Rat Serum	30	3.39-149.3 3	100	31.78
Monkey Serum	2	9.74-10.17	100	9.95

n.d. = non-detectable. Samples measured below the sensitivity are considered to be non-detectable.