

# Mouse VEGF ELISA Instructions

## Cat: EM0013

### Content

	CAT	Volume
1 CP (Coated Plate)	EM0013CP	96 well
2 S (Standard)	EM0013S,S1~S7,S0	9 vial
3 DA (Detect Antibody)	EM0013DA	6 ml/bottle
4 SH (Streptavidin-HRP)	ESH01	12 ml/bottle
5 AB (Assay Buffer 1×)	EAB01	12 ml/bottle
6 TS (TMB Substrate)	ETS01	12 ml/bottle
7 SS (Stop Solution)	ESS01	12 ml/bottle
8 WB (Wash Buffer 10×)	EWB01	50 ml/bottle
9 SF (Sealer Film)	ESF01	6 piece

**NOTE:** After the kit is opened, the stabilization period of each content is 30 days, so please use the kit within 30 days after opening.

### REAGENT PREPARATION

#### Washing Buffer (1×) Preparation


Pour entire contents (50 ml) of the **Washing Buffer Concentrate** (10×) into a clean 500 ml graduated cylinder. Bring to final volume of 500 ml with glass-distilled or deionized water. Transfer to a clean wash bottle and store at 2 to 25°C.

#### Standard Curve Preparation:

S1 to S7 and S0 is ready to use for serum and plasma.

Other sample type, prepare the standard curve with whatever buffer (SPB, Sample Prepared Buffer) is used to prepare the sample, such as cell culture supernatant, tissue grinding liquid, cell lysate, etc. Urine sample use AB (Assay Buffer) prepare standard curve.

The Mouse VEGF Standard EM0013S 30000 pg/ml 30  $\mu$ l + 270  $\mu$ l SPB serves as the high standard (3000 pg/ml). Pipette 200  $\mu$ l of SPB into each tube. Use the high standard to produce a 1:2 dilution series. Mix each tube thoroughly before the next transfer. SPB serves as the zero standard (0 pg/ml).

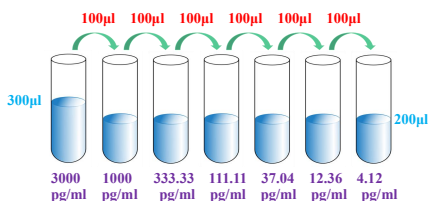


### ASSAY PROCEDURE

Bring all reagents and samples to room temperature before use.

- 1 Prepare all reagents and working standards as directed in the previous sections.
- 2 Remove excess **CP** (Coated Plate) strips from the plate frame, return them to the foil pouch and reseal.
- 3 Add 50  $\mu$ l of **AB** (Assay Buffer) to each well.
- 4 Add 10  $\mu$ l of **Standard or sample** per well. Ensure reagent addition is uninterrupted and completed within 15 minutes.
- 5 Add 50  $\mu$ l of **DA** (Detect Antibody) to each well.
- 6 Cover with an **SF** (Sealer Film). Incubate at room temperature (18 to 25°C) for 1 hours on a microplate **shaker** set at 500 rpm.
- 7 Aspirate each well and **wash**, repeating the process four times. Wash by filling each well with WB (Washing Buffer 300  $\mu$ l). Complete removal of liquid at each step is essential to good performance. After the last wash, remove any remaining WB (Washing Buffer) by aspirating or decanting. Invert the plate and **blot** it against clean paper towels.
- 8 Add 100  $\mu$ l of **SH** (Streptavidin-HRP) to each well.
- 9 Cover with a new **SF** (Sealer Film). Incubate at room temperature (18 to 25°C) for 30 min on a microplate **shaker** set at 500 rpm.
- 10 Repeat aspiration/**wash** as in step 7.
- 11 Add 100  $\mu$ l of **TS** (TMB Substrate) to each well. Incubate for 5-30 minutes at room temperature.
- 12 Add 100  $\mu$ l of **SS** (Stop Solution) to each well.
- 13 Determine the optical density within 30 minutes, using microplate **reader** set to 450 nm corrected with 570 nm or 630 nm.

## TYPICAL DATA

Mouse VEGF Typical Standard

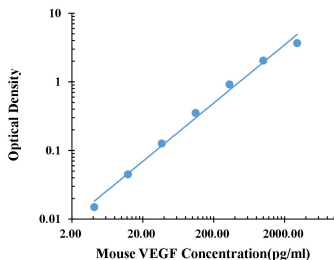

pg/ml	O.D.	Average	Corrected
0.00	0.0084	0.0091	0.0088
4.12	0.0241	0.0231	0.0236
12.35	0.0557	0.0526	0.0542
37.04	0.1373	0.1337	0.1355
111.11	0.3702	0.3549	0.3626
333.33	0.9362	0.9147	0.9255
1000.00	2.0741	2.0273	2.0507
3000.00	3.6762	3.6825	3.6794

## SENSITIVITY

The minimum detectable dose (MDD) of mouse VEGF is typically less than 0.57 pg/ml.

The MDD was determined by adding two standard deviations to the mean optical density value of ten zero standard replicates and calculating the corresponding concentration.

## PRECISION

**Intra-assay Precision (Precision within an assay)** Three samples of known concentration were tested twenty times on one plate to assess intra-assay precision.

**Inter-assay Precision (Precision between assays)**

Sample Number	Intra-assay Precision			Inter-assay Precision		
	S1	S2	S3	S1	S2	S3
	22	22	22	6	6	6
Average (pg/ml)	59.4	340.6	878.4	64.3	347.1	993.3
Standard deviation	1.8	9.8	21.3	2.7	12.1	20.0
Coefficient of variation (%)	3.0	2.9	2.2	4.2	3.5	2.0

## RECOVERY

The spike recovery was evaluated by spiking 3 levels of mouse VEGF into health mouse serum sample. The un-spiked serum was used as blank in this experiment.

The recovery ranged from 114% to 131% with an overall mean recovery of 122%.

## LINEARITY

To assess the linearity of the assay, five samples were spiked with high concentration of VEGF in mouse serum and diluted with Sample Diluent to produce samples with values within the dynamic range of the assay.

The linearity ranged from 106% to 124% with an overall mean recovery of 115%.

## SAMPLE VALUES

Serum/Plasma – Thirty samples from apparently healthy mice were evaluated for the presence of VEGF in this assay. No medical histories were available for the donors.

Sample Matrix	Sample Evaluated	Range (pg/ml)	Detectable %	Mean of Detectable (pg/ml)
Serum	30	16.98-45.87	100	27.58

n.d. = non-detectable. Samples measured below the sensitivity are considered to be non-detectable.